

# eyeForce

Item No.: 1800 0000

Our most innovative and modern patient positioning system with auto-run functions - specialized for precision surgery in head area. With optional *e-drive* easy-push aid, Steering Guide and SES automated extension system - for even more ergonomics and economy at your workplace.

# Trendsetting comfort and supreme quality



# MED.TEC at its highest level

#### **Double castors**

Four fully rotating double-run castors for back-friendly and easy rolling motion.

#### Chassis

Precision lightweight chassis made of cast aluminum.

#### Frame

Lasercut procedure. Highly stable stainless steel frames. Including standard rails 25 x 10 mm for the attachment of a wide range of components.







14

# Precise positioning of the head part

## Sterile adjustment by surgeon/maximum legroom

#### Headrest

3D and longitudinally adjustable headrest element with head molding and neck support for guaranteed safe and correct positioning of the head.








## Components

#### Headrest with ball joint and more

Discover additional features and extensive components on page 22-23.

Item No.: 3351 1020

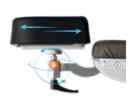





### Tiltable foot control

Tiltable foot keypad via two 4-way joysticks plus activation switch allow sensitive and precise chair adjustment even after the surgeon is already steril. Pushing the pedal lip (by foot) will activate the electronic control and battery mode of the equipment.



()

# Positioning

The *eyeForce* treatment chair is adapted for patient positioning and ergonomic and comfortable patient treatment, in particular for outpatient surgical techniques, and offers versatile adjustment of the head, back, seat and leg sections and overall chair height for patient transport.

• Five freely configurable memory positions

Positioning by hand control

Standard recommended memo keys

- Reset button
- Auto run function
- Battery charge indicator



## Hand control

#### Activate memory function/automatic movement

It is possible to save 5 MEMO positions: Set the desired position (hold down the MEMO key until the orange LED lights up). Then press the numeric key of a memory location for a MEMO position - the MEMO position is now saved.

Press and hold "AUTO" until the LED lights orange, then briefly tap on either 1, 2, 3, 4 or 5 to move to the desired position of the surgical bed.

AUTO

All castors are braked (red)



# Transport Anesthesia

# "Steering Guide" maneuvering aid

#### Operation via push handles and hand control

The *eyeForce* treatment chair offers an electronically controlled braking system with "Steering Guide" function. In addition to "free movement" and "brake applied" the treatment chair offers directional guidance as an additional electronic wheel positioning function. This additional function facilitates the controlled maneuvering of the chair by just one person and can be activated from both sides of the chair via the chair handles or optionally via the hand control.

Directional guidance is activated (orange)

Directional guidance is deactivated (green)



# Surface-Extensions-System "SES" positioning aid

#### Automated solution for patient positioning

Both the backrest and the legrest segments automatically move 50 mm along with the patient when moving the chair into the reclined or entry position. The "SES" avoids patient repositioning, which would otherwise be necessary, by extending and flexing the individual segments. Returning the leg segment to the entry position reduces the entry height by an additional 50 mm.


#### **Retractable push handles**

Manually adjustable push handles on the foot end with automatic return function, for easy patient transfer.



#### Movement in reclining position

Pushing from the head direction

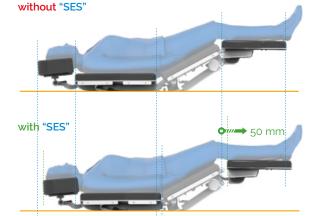

**0** *m* → 50 mm



#### without "SES"

Movement in entry position

with "SES"



# eyeForce + e-drive

Item No.: 1800 0000 + 3316 0000

#### Powerful yet gentle.

()

EMIUM

*e-drive* system for effortless and rapid movement of the chair over long distances. The new optional system consists of a 5th wheel with drive motor. It serves as a pushing aid for easy forward and backward movement – with an extra range of 2 kilometers per battery charge. The system's operation is identical to that of the "Steering Guide" system, via the respective push handles with potentiometer for infinite and precise speed adjustment. It is possible to integrate the system any time at a later date.



## The electronic push-aid

#### Control over push handles

Push handle at the left head/foot end with a rocker switch for moving the treatment chair forward and back.

Push handle at the right head/foot end with a potentiometer for speed control.






## Low-noise movement/trackless movement

By switching on the fifth electronically drive wheel in the chassis of the chair, the physical labor of the person doing the pushing is noticeably reduced.

- Shock-free transport
- Long-lived design
- Optimum floor contact
- Long range

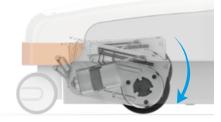

# The ease of pushing - How much force is needed to move 300 kg? - Zero.

#### Who would have thought this is possible?

This is how far your staff pushes


- a load of 300 kilograms per year:
- approximately 70 meters per surgery round
- up to 10,000 surgeries per year each

This results in up to **700 km per year!** 



# Push handles 2.0 - perfect interaction - top-of-the line mobility and ergonomics.



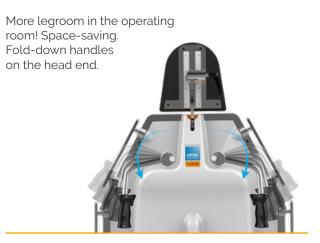



# **Technical specifications**

#### Chassis

- Outside length
- Outside width
- $\cdot\,$  Brake system acting on all four castors
- Directionally guided steering ("Steering Guide"), easily activated at the push of a button
- Castor diameter, electr. conductive ø 150/100 mm
- Foot control element with joystick (total, head, back, Trendelenburg up/down)
- Emergency stop button
- Optional push-aid (e-drive)
- With equipotential bonding connection
- Abrasion-resistant corner guards (orange) made of synthetic rubber

| Upper frame |
|-------------|
|-------------|

900 mm

662 mm

| <ul> <li>4-part support frame (head, back, seat, leg part)</li> </ul> |               |                    |  |  |  |  |
|---|---------------|--------------------|--|--|--|--|
| Removable uphols  | 40 mm         |                    |  |  |  |  |
| • Width of upholste   | 650 mm        |                    |  |  |  |  |
| Standard stainless  | steel side ra | ils                |  |  |  |  |
| on both sides   | 25 x 10 mm    |                    |  |  |  |  |
| • Width incl. norm ra   | 700 mm        |                    |  |  |  |  |
| • Length  |               | 1,780 – 1,900 mm   |  |  |  |  |
|   |               |                    |  |  |  |  |
|   | Angle         | Dimensions (L x W) |  |  |  |  |
| <ul> <li>Head part</li> </ul>   | -31°/+22°     | 255 x 245 mm       |  |  |  |  |
| <ul> <li>Back part</li> </ul>   | 0°/+75°       | 540 x 650 mm       |  |  |  |  |
| <ul> <li>Seat part</li> </ul>   | -6°⁄+26°      | 450 x 650 mm       |  |  |  |  |
| <ul> <li>Leg part</li> </ul>  | 0°⁄-78°       | 510 x 610 mm       |  |  |  |  |
| <ul> <li>Height adjustment</li> </ul>                                 | t             | 640–900 mm         |  |  |  |  |
| <ul> <li>Entry height</li> </ul>                                      |               | 580 mm             |  |  |  |  |
| Legroom   |               | 470 – 730 mm       |  |  |  |  |
| Trendelenburg -26°/Reverse Trendelenburg +6°                          |               |                    |  |  |  |  |

#### Electrical control

#### Manual keypad

- Memory controller with 10 buttons for head, back, seat, leg part – total up/down; 2 buttons for directionally guided steering, 1 button for electronic brake, 1 easy entry button, memory button, on-off button, auto-run function
- 5 freely programmable memory positions
- Visual battery level indicator

#### Foot control element

- Foot control element 2 joysticks plus activation switch
- Optional push-aid (e-drive)

#### Electric drives

- Head, back, seat, leg part, total up/down
- Trendelenburg/Reverse Trendelenburg



#### Powerful accumulators

- Power supply-independent battery unit with 2 batteries (lead-gel), pollution free, maintenance-free, cadmium-free
- Quick-exchange battery drawer replaceable without tools

| <ul> <li>Universal battery charger</li> </ul> | Input:  | 100–240 V |
|---|---------|-----------|
|   | Output: | 24 V      |

Operating time per battery charge

| <ul> <li>Measured by operating circle<br/>(Prepare-Operate-Recover)</li> <li>Quantity of operating circles depends<br/>on individual operation</li> </ul> | 100 OPs |
|---|---------|
| Max. load   | 250 kg  |
| Weight  | 163 kg  |
|   |         |

Subject to technical modifications.

#### Mains-free operation

Easy recharging of the equipment thanks to charging adapter and integrated battery charger circuit board.

#### Simple battery replacement

An additional quick-exchange drawer (including batteries) allows you to quickly replace the rechargeable battery packs (back-up system **Item No.**: *3339 0000*).





#### **Equipment protection**

#### Cover drape stayClean

2,000 x 730 mm, box á 160 Pieces 4x 40 Pieces **Item No.:** *3034 0014* 

Page 64





#### Upholstery/color card

Standard color

anthracite Item No.: 3037 0052

#### **"skai Tundra"** Page 102

**"skai Pandoria Plus"** Page 103



# **Components**



•



eyeForce Item No.: 1800 0000



e-drive, upcharge

Item No.: 3316 0000



Backup battery drawer incl. batteries

Item No.: 3339 0000

Mounting clamp polished

incl. locking device

Item No.: 3090 0000

aluminium, with guick clamp

lever for norm rail 25 x 10 mm,



Headrest with ball joint, 360° all-round adjustment, factory-fitted

Item No.: 3351 1020

Flexible surgical drape

3090 0000 not included)

Item No.: 3032 0000

support rod (mounting clamp



incl. guick-release clamp for

mounting on norm rail

Item No.: 3351 0050

580 x 200 mm,



Side rails right and left, swiveling incl. locking, upholstered



Hand control holder, PVC (mounting clamp 3090 0000 not included)

Item No.: 3018 0009



Anesthesia arm board

with universal joint, incl. mounting clamp

Item No.: 3011 0000


Hand control holder for attachment onto an I.V. Pole Ø 25/18 mm, PVC

Item No.: 3027 0006



Head part bumper stainless steel

Item No.: 3028 0000

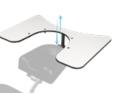

Wrist rest 550 x 370 mm, height adjustable, incl. mounting plate, adjustment 115 mm



Flexible surgical drape

connection (mounting clamp

3090 0000 not included)

support rod, with air

Item No.: 3033 0000

Wrist rest, adjustable, incl. mounting plate 3013 0000, adjustment 115 mm

Item No.: 3014 0000



Wrist rest, adjustable, tiltable, incl. mounting plate 3013 0000, adjustment 115 mm

Item No.: 3014 0005

|   |  |   |  |   | Hygie  | ne products  |
|---|--|---|--|---|--|--|
|   |  |   |  |   | Å  |  |
| Mounting plate stainless steel,<br>for wrist rest 3014 0000,<br>3014 0005 and 3209 0000 | Head fixation pillow with<br>hollow (rubins pillow),<br>410x250 mm, height 100 mm,<br>factory-fitted | Shoulder parts,<br>right and left                           | Instrument table<br>500×300×15 mm,<br>height adjustable, swiveling,<br>detachable,<br>ingl. according to the | Component arm<br>stainless steel, swiveling,<br>incl. extendible I.V. Pole<br>Ø 25 mm | Cover drape for  | Kick bucket,   |
| Item No.: 3013 0000   | Item No.: 3022 0000  | Item No.: 3043 0000   | incl. component arm<br>Item No.: 3065 0003   | Item No.: 3024 0000   | treatment chair,<br>stay <b>Clean</b><br>2,000 x 730 mm,<br>box á 160 pieces | capacity 12l,<br>stay <b>Clean</b> kick bucket                                     |
|   |  | AG  |  |   | ltem No.: 3034 0014  | ltem No.: 3351 1046  |
| Wrist secure straps,<br>pair  | Body strap silicone<br>incl. clamp, 3-piece  | Foot protection foil  | Foot part extension,<br>attachable, incl. norm rail  | Step  |  |  |
| Item No.: 3015 0000   | ltem No.: 3351 1044  | ltem No.: 3328 0000   | Item No.: 3351 0002  | Item No.: 3351 0042   | Medic  | al product   |
|   |  |   |  |   |  | -1-  |
| Knee half roll,<br>600x280x140mm,<br>PUD-coated   | Knee roll round,<br>500 x 150 mm,<br>PUD-coated  | Upholstery incl. integrated<br>knee roll,<br>factory-fitted | Cylinder holder  |   |  | Instrument table<br><i>support<b>Line</b></i><br>height adjustable<br>730–1,000 mm |

Item No.: 3072 0001

Item No.: 3351 0026

Item No.: 3050 0000

Item No.: 3019 0000

Item No.: 3037 0266

**PREMIUM** LINE  $\bigcirc$ 

# "skai<sup>®</sup>Tundra"

High quality skai® material with fine calf leather grain and subdued matte-glosseffect. The very extensive and trendorientated color spectrum enables many applications and combinations.

| <ul> <li>Product Specification</li> <li>Long-lasting</li> <li>Tear resistant</li> </ul>   |   | weiß F6461198<br>Item No.:134    | creme F6461210<br>Item No.:057         | beige F6461636<br>Item No.:094          | sand F6461230<br>Item No.:058            | pistazio F6461610<br>Item No.: <i>153</i> | mint F6461495<br>Item No.:067        | aqua F6461456<br>Item No.:001   | türkis F6461338<br>Item No.:008 |
|---|---|----------------------------------|--|---|--|---|--------------------------------------|---|---------------------------------|
| <ul> <li>Easy to clean</li> <li>Abrasion resistant (DIN EN</li> <li>Lightfast</li> <li>Disinfectant resistant</li> <li>Biocompatibility:<br/>cytotoxicity (EN ISO 10 993</li> </ul> | 3-5) and  | perle F6461288<br>Item No.:007   | silber F6461550<br>Item No.:070        | stone F6461608<br>Item No.:075          | bleu F6461496<br>Item No.:068            | dolphin F6461609<br>Item No.:076          | sky F6461581<br>Item No.: <i>073</i> | royal F6461196<br>Item No.:051  | atoll F6461453<br>Item No.:062  |
| epicutaneous test (EN ISO<br>• Blood- and urine resistant<br>• Phthalate-free<br>• Vegan  |   | violett F6461356                 | cyclame F6461381                       | orchidee F6461354                       | bordeaux F6461350                        | kirsche F6461199                          | feuer F6461454                       | karneol F6461557  | orange F6461556                 |
| eye <b>Force</b><br>500 <b>XLE</b><br>500 <b>XLE</b> comfort  | 1800 0000<br>1500 0006<br>1500 0006<br>† 3037 0080            | Item No.:059                     | Item No.:119                           | Item No.: 011                           | Item No.:010                             | Item No.:054                              | Item No.:063                         | Item No.:072  | Item No.:071                    |
| 500 XLS<br>500 ECO<br>400 XL<br>surgiForce Formed seat<br>Formed seat   | 1500 0005<br>1500 0008<br>1400 0001<br>2600 0008<br>2600 0009 | apricot F6461458<br>Item No.:066 | mais F6461457<br>Item No.:065          | mango F6461607<br>Item No.:074          | safran F6461542<br>Item No.:069          | limonade F6461343<br>Item No.:009         | limone F6461543<br>Item No.:109      | smaragd F6461455<br>Item No.:064  | agave F6461273<br>Item No.:006  |
| surgiLine Formed seat<br>assistLine<br>motionLine<br>motionLine-E   | 2200 0008<br>2400 0015<br>1700 0000<br>1700 0001              | chrom F6461205<br>Item No.:056   | ocean F6461204<br>Item No:: <i>055</i> | schoko F6461357<br>Item No.: <i>060</i> | schwarz F6461195<br>Item No.: <i>050</i> | anthrazit F6461197<br>Item No.:052        |                                      | The colors displaye<br>may differ slightly fr<br>colors. If you wish to | rom the original                |

+ Item No.: 3037 0...

colors. If you wish to receive original color charts, please contact us at info@ufsk-osys.com.





# "skai<sup>®</sup> Pandoria Plus"

skai<sup>®</sup> Pandoria Plus is a high-quality upholstery material, whose fine micropore structure gives a certain unobtrusiveness, which is refined by an attractive tone-in-tone effect. The very pleasant feel of the material and the excellent material characteristic and the broad, modern color spectrum make skai<sup>®</sup> Pandoria Plus universally usable.

#### **Product Specification**

Long-lasting

()

- Tear resistant
- Easy to clean
- Abrasion resistant (DIN EN ISO 5470-2)
- Resistant to oil and fat
- Lightfast
- Disinfectant resistant
- Flame protection
- Biocompatibility: cytotoxicity (EN ISO 10 993-5) and epicutaneous test (EN ISO 10 993-10)
- Resistant to sudor, blood and urine
- Antibacterial according to JIS Z 2801
- Vegan

| eye <b>Force</b>         |             | 1800 0000   |
|--------------------------|-------------|-------------|
| 500 <b>XLE</b>           |             | 1500 0006   |
| 500 <b>XLE</b> comfort   |             | 1500 0006   |
|                          |             | + 3037 0080 |
| 500 <b>XLS</b>           |             | 1500 0005   |
| 500 <b>ECO</b>           |             | 1500 0008   |
| 400 <b>XL</b>            |             | 1400 0001   |
| surgi <b>Force</b>       | Formed seat | 2600 0008   |
| surgi <b>Force</b> light | Formed seat | 2600 0009   |
| surgi <b>Line</b>        | Formed seat | 2200 0008   |
| assist <b>Line</b>       |             | 2400 0015   |
| motion <b>Line</b>       |             | 1700 0000   |
| motion <b>Line</b> -E    |             | 1700 0001   |
|                          |             |             |

+ Item No.: 3037 0136

+ Item No.: 3037 0...



The colors displayed on your device may differ slightly from the original colors. If you wish to receive original color charts, please contact us at info@ufsk-osys.com.

# Upgrade your OR to perfection.

